

Quad-Tray® XL

Xtra Long



The **Quad-Tray® XL (Xtra Long)** incorporates all the benefits of the original **Quad-Tray** from **Clinician's Choice®**, which includes an inflexible aluminum construction, wide, customizable arch, and thin distal bar. However, the **Quad-Tray XL** was designed with two additional benefits in mind; a full quadrant design, to assist the dental lab in providing ideal occlusion by recording cuspid rise, and more lateral support for the impression material, increasing rigidity and reducing impression flex.

Full Quadrant Design

The **Quad-Tray XL** allows dentists to take ideal, distortion-free dual-arch impressions in cases that involve multiple preparations, including small bridges. Both the lingual and buccal walls of the **Quad-Tray XL** are enhanced to provide a solid lock of the impression material to the tray.


On the buccal side, as the tray curves towards the anterior, the **Quad-Tray XL's** perforated walls provide sufficient height to encompass the height of anterior teeth.






Improved Occlusion

The more occlusal information the laboratory receives from your impression, the more precise the occlusion that can be generated. Whether you are impressing one or multiple teeth, the **Quad-Tray XL** will provide more detail which is a significant clinical advantage.

Extended Length

The extended length also allows the clinician to capture larger cases, including 3-unit bridges, while providing greater lateral support for the impression material, thereby reducing the risk of distortion.

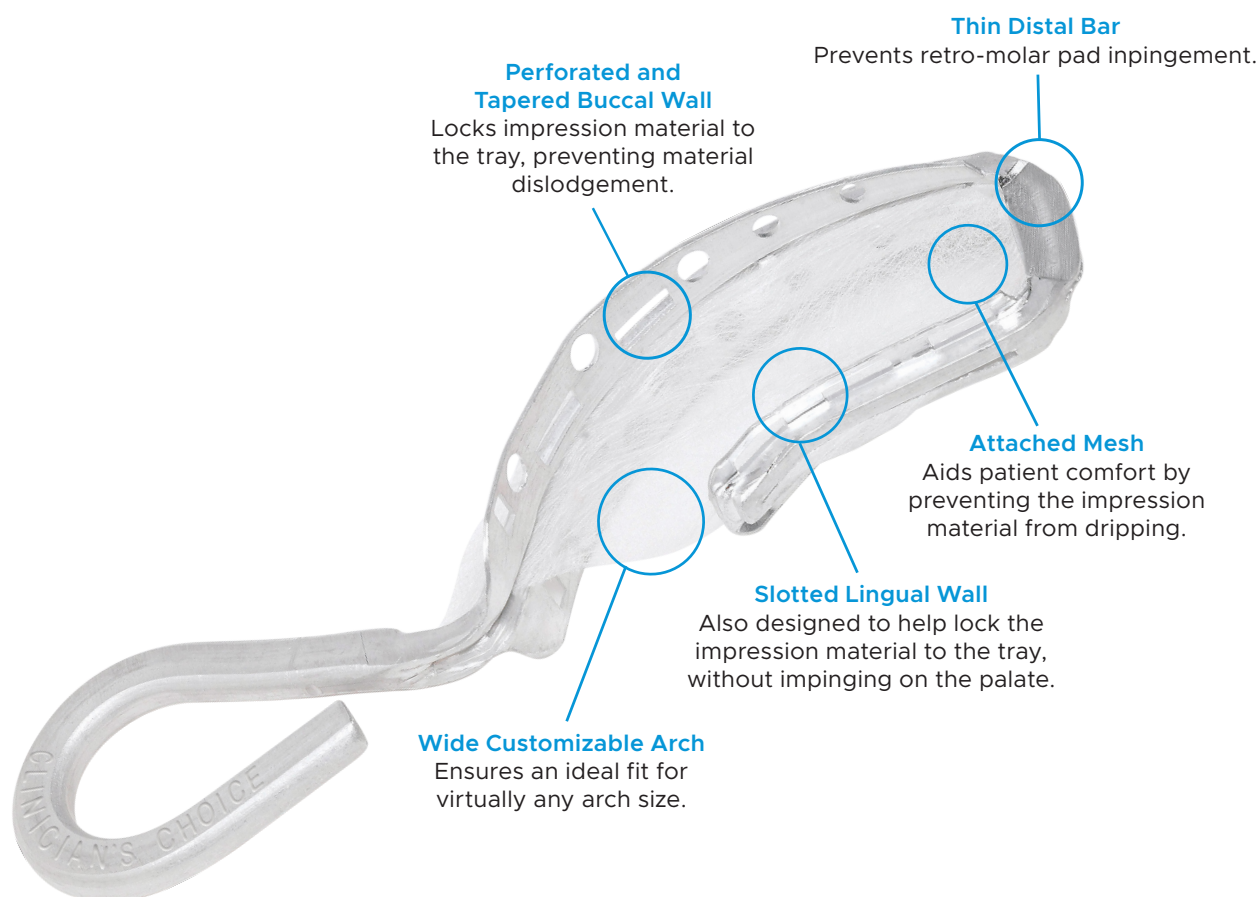

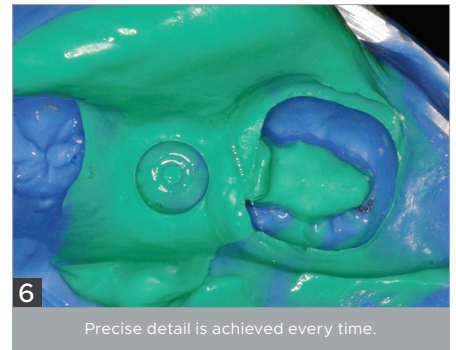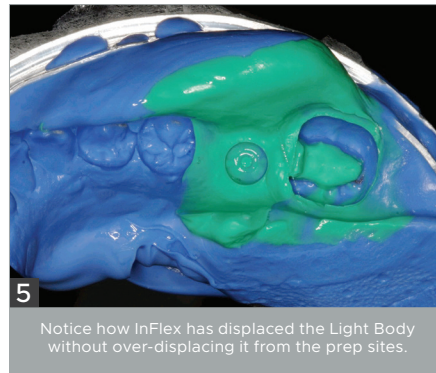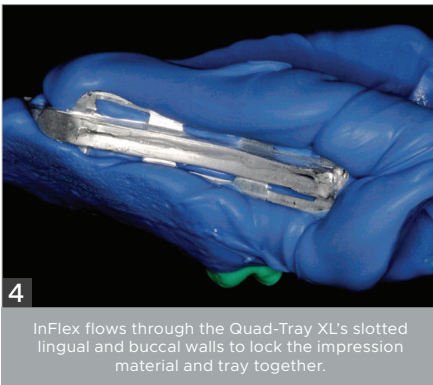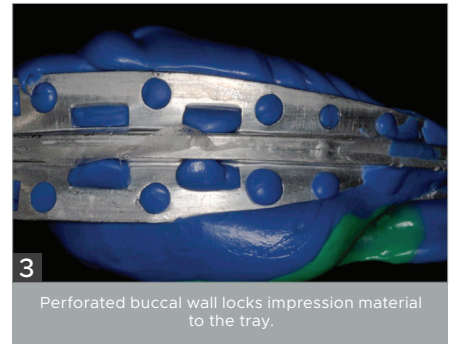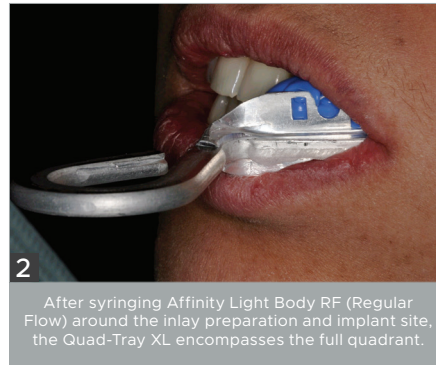

QUAD-TRAY XL 40-PACK (041400)

QUAD-TRAY TRIAL KIT (040152)

Contains: 20 x Quad-Tray Xtreme,
10 x Anterior Quad-Tray X2,
20 x Quad-Tray XL